



TOP
INNOVATION
2017 - 2018

Clotting
Chromogenic
Immunturbidimetric

Coatron

Semi-automated
Coagulation Analyzer Series

With 1, 2 or 4 optical channels



TECO
Innovation in Coagulation

A new area of manual and semi-automated Coagulation Analyser rise up

The Coatron X instrument line is a consequent continuation in the development of the Coatron product line. Over 25 years in experience and innovation is the reference for our new Coatron X instrumentation line.

The unique detection principle in combination with the high-level analytical algorithm calculates exact, precise and reproducible results.

Easy in operation – self instructing user dialogue - reliable

Highest optical resolution, enlarged optic range, smallest sample and reagent volume

0.1 mOD, 0 - 3800 mOD, just with 75 µL sample and reagent volume

Complete optical analysis

No further parts required, like balls, stirrers etc.

Adaptation of the light level

Automatic light level adjustment of the optic channels to each sample

Exclusion of disturbance

Stray light reduction, exact temperature control, all parameter are preset

“Complete range of Coagulation Analysis with the highest standard and reliability. The new generation of Coagulation instruments with optical detection are here.”

Coatron X - product family

With 1, 2 or 4 optical channels.

www.teco-medical.com

Prepared for the daily routine and the upcoming requirements

One instrument – many possibilities

The Coatron X family is prepared to work with one, two or four channels. The built-up and functionality is specifically designed to each instrument version and requirements. The operation with the intuitive user dialogue and handling of the detection results are easy and effective. The possibility to connect the instrument to the **TECO Cloud** offers new perspective of instrument, reagent and consumables verification and handling. The precise and correct patient result is what we want to secure.



Quality is our basic demand

TECO develop and produce with qualified and specialized companies, located in Germany. High reliability, nearly maintenance free instruments are our benefit. Our reference is 25 years, in worldwide laboratories, with satisfied users.



TECO Cloud Services – A strong data bank and application service behind

All instrument versions of the Coatron X family are connectable via Bluetooth to Smart-devices, like mobile devices, tablets, etc. with a specific APP or direct access to the TECO Cloud Services.



Coatron	Eco	Pro	Top
General			
Dimensions	230 x 148 x 94 mm (l, b, h)		
Display	Colored touch display 4.3"		
Pre-warm temperature	37°C		
Pre-warm cuvettes (pcs.)	10	20	20
Pre-warm reagent 24mm (pcs.)	1	1	1
Pre-warm reagent 22mm (pcs.)	2	2	2
Pre-warm reagent 11mm (pcs.)	2	2	2
Reagent mixing position	-	1	1
Power values	110-240Vac, 50-60Hz / 5Vdc, 3.3A		
Interfaces			
RS232 (2x)	Printer, Barcode reader		
USB (2x)	Network, Firmware update		
Bluetooth	TECO Cloud, App		
Optic / tests			
Optic channels	1	2	4
Wavelength (nm)	620 (red)	405 (UV)	405 (UV)
Global Coag. tests	PT, APTT, TT, FIB		
Specific Coag. tests	-	individual factors	
Chromogenic Coag. tests	-	AT, PC	
Latex based tests	D-Dimer		
Whole blood tests	PT-B	-	



The details make the difference

Coatron X

The remarkable details in every single component is achieved by selecting of premium suppliers.

The performance of a high level instrument is strongly depending on the concept in general and the perfect usability to reach the requirements of a modern laboratory analyser.

Priority No. 1 was to get a daily routine reliability and easy-to-use operation.

Software and connection possibilities

With the Coatron X product family starts a new time line in analysis management and service maintenance. Operation via intuitive, colored touchscreen, as well patient result management are perfectly optimized.



Operation details

Coatron	Eco	Pro	Top
Operation			
Touchscreen 4.3"	✓	✓	✓
Real time clock	✓	✓	✓
Stopwatch	✓	✓	✓
Language selection	✓	✓	✓
Interfaces			
USB to LIS	✓	✓	✓
Network to LIS (TECAM software required)	✓	✓	✓
Management			
Test calibration	✓	✓	✓
Tracking to Pat.ID, Patient ID, Sample ID or Auto ID	✓	✓	✓
Automatic optic start (no Starterpipette required)	✓	✓	✓
Double determination	✗	✓	✓
Sample management (ID)	✗	✓	✓
Reagent management (ID) (lot und expiry)	✗	✓	✓
Internal result databank	✗	✓	✓
Patient identification with barcode		optional	



Intuitive operation and control

Clear and easy to operate user dialogue with a high quality colored touchscreen

- Direct usable
- Short learning phase
- Logic, intuitive operation
- Reliable touchscreen surface
- Quick touch response



For small and mediate laboratory requirements

Concept is suitable for daily routine work in Coagulation laboratories and hospitals

- Three different versions available, depending on number of samples per day
- In maximum up to 4 optic channels available

Interfaces

RS232 (2x)

- For external serial printer and external barcodereader

LIS/USB

Bluetooth



Integrated barcode scan for reagents.






TECO Cloud Services

A strong data bank and application service behind

All instrument versions of the Coatron X family are connectable via Bluetooth to Smart-devices, like mobile devices, tablets, etc. with a specific APP or direct access to the TECO Cloud Services.

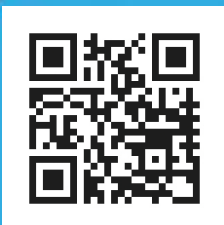

For trading partners worldwide, please visit our web-page

TECO Medical Instruments Production + Trading GmbH
Dieselstr. 1, 84088 Neufahrn, Germany
Tel.: +49 (0) 8773 70780-0, Fax +49 (0) 8773 70780-29
info@teco-gmbh.com, www.teco-medical.com

TECO
Innovation in Coagulation